

KV β 1 Rabbit pAb Catalog NO.: EA307 For research use only.

Overview

Product name KV β 1 Rabbit polyclonal antibody

Source Rabbit

Applications IHC

Species reactivity Human, Rat, Mouse

Recommended dilutions Immunohistochemistry:1/100-200

NOTE: Optimal dilutions should be determined by the end user.

Immunogen Synthetic Peptide

Species Human

Storage PBS with 0.02% sodium azide and 50% glycerol pH 7.4.

Store at -20° C. Avoid repeated freeze-thaw cycles.

Isotype IgG

Clonality Polyclonal

Concentration 1 mg/ml

Observed band 50kDa


GenelD (Human) 7881

Human Swiss-Prot No. Q14722

Cellular localization Cell membrane, Cytoplasm, Membrane

Alternative Names hKvb3,KCNA1B,KV beta 1,Voltage-gated potassium channel beta 1 subunit

Background N/A

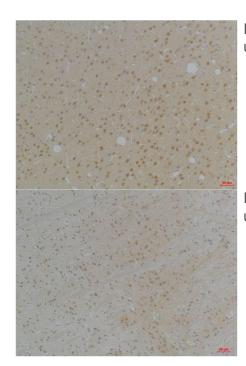

Immunohistochemical analysis of paraffin-embedded Rat BrainTissue using Kvb1 (EA307) Rabbit pAb diluted at 1:200.

Immunohistochemical analysis of paraffin-embedded Mouse BrainTissue using Kvb1 (EA307) Rabbit pAb diluted at 1:200.